


YOUR ALTIUS CVC RANGE



REDUCE
PATIENT RISKS
WITH ALTIUS

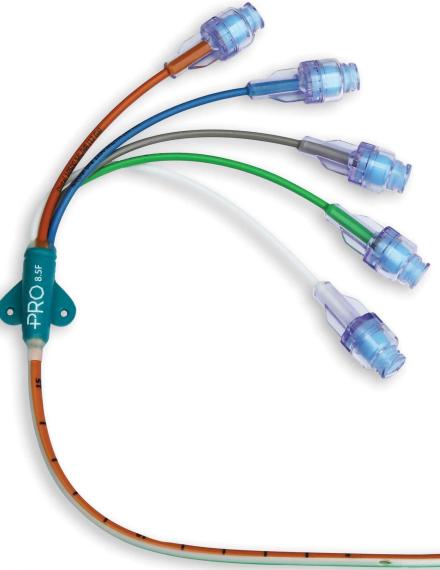

WHY CHOOSE **ALTIUS™?**



Altius ProActiv+™ is a chlorhexidine free.

non-eluting, covalently bonded antimicrobial to

protect your patient for the duration of use.

patient's bloodstream and sustained antimicrobial activity.

Altius® is designed based on the feedback received from healthcare professionals for 21st century healthcare.

Altius is manufactured utilising our proprietary Multi-Tube™ technology to achieve perfectly round lumens and laminar flow.

THE RANGE

• Adults: 1 to 5 and 7 lumen.

For French size, length and flow rates please see ordering guide.

















5ml Guidewire Introducer Luer Slip Syringe: Facilitates



Altius Guidewire Dispenser: Ergonomic design to facilitate guidewire advancement even with wet gloves.



Kink-Resistant Nitinol Guidewire: To prevent guidewire kinking even in difficult anatomy.




No. 11 Mini Scalpel: To make a small incision prior to dilation.



Hydrophilic Dilator: For smooth dilation to prevent vessel trauma.



Secondary **Fixation Wing**

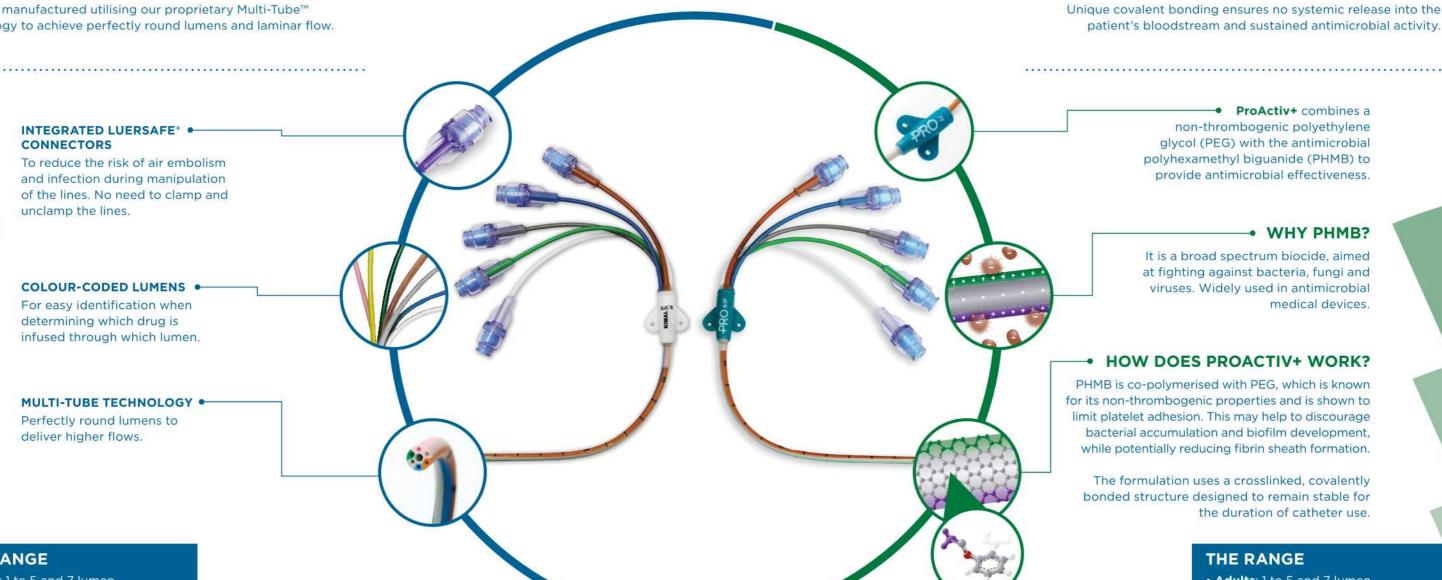

 ProActiv+ combines a non-thrombogenic polyethylene glycol (PEG) with the antimicrobial polyhexamethyl biguanide (PHMB) to provide antimicrobial effectiveness.

WHY PHMB?

It is a broad spectrum biocide, aimed at fighting against bacteria, fungi and viruses. Widely used in antimicrobial medical devices.

HOW DOES PROACTIV+ WORK?

PHMB is co-polymerised with PEG, which is known for its non-thrombogenic properties and is shown to limit platelet adhesion. This may help to discourage bacterial accumulation and biofilm development, while potentially reducing fibrin sheath formation.

The formulation uses a crosslinked, covalently bonded structure designed to remain stable for the duration of catheter use.

THE RANGE

· Adults: 1 to 5 and 7 lumen. For French size, length and flow rates

please see ordering guide.













The insertion kit has been designed to reduce the risks during insertion.








guidewire advancement.



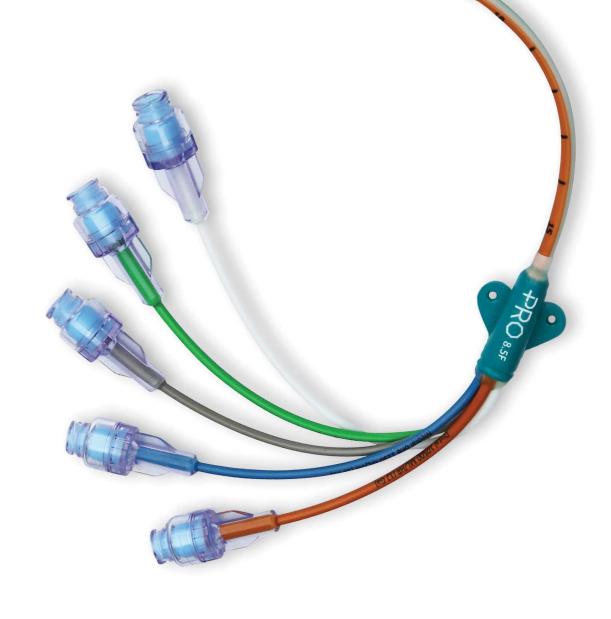

For more information, pricing or samples, speak to your local Kimal Representative, or contact:

UK

Tel: +44 (0)845 437 9542 Email: customercare@kimal.com INTERNATIONAL

Tel: +44 (0)845 437 9543 Email: internationalenquiries@kimal.com **MIDDLE EAST, ASIA & AFRICA**

Tel: +97 16 552 8146 Email: mea@kimal.com

© 2025 Kimal PLC. All rights reserved. No part of this document may be reproduced in any form or by any means without permission in writing from Kimal PLC.

Kimal PLC reserve the right to change the specification of its products at any time.

Kimal PLC, Clayfield Rd, Worcester WR4 0AE.

